



Ordering Information (cont.)

Cat. No.	Description	Size	Each/Case	
3M™ Synthetic Cast Stockinet				
MS01		1" x 25 yds. (2,5 cm x 22,8 m)	1	
MS02		2" x 25 yds. (5 cm x 22,8 m)	1	
MS03		3" x 25 yds. (7,5 cm x 22,8 m)	1	
MS04		4" x 25 yds. (10,1 cm x 22,8 m)	1	
MS06		6" x 25 yds. (15,2 cm x 22,8 m)	1	
MS08		8" x 25 yds. (20,3 cm x 22,8 m)	1	
3M™ Synthetic Cast Padding			Each/Bag	Bag/Case
CMW02		2" x 4 yds. (5 cm x 3,6 m)	20	4
CMW03		3" x 4 yds. (7,5 cm x 3,6 m)	20	4
CMW04		4" x 4 yds. (10 cm x 3,6 m)	20	4
CMW06		6" x 4 yds. (15 cm x 3,6 m)	20	4

For additional information contact your 3M Health Care Sales Representative, call the 3M Health Care Customer Helpline at 1-800-228-3957, or visit our web site at www.3m.com/immobilization.

Ordering Information

IMPORTANT: When ordering 3M™ Scotchcast™ products, please ensure that you quote the correct code number. Your chosen roll width must replace the asterisk (*) when ordering. For example, 8200*L would become 82003L for 3" Light Blue material and 4" Red material would be 82004R.

Cat. No.	Description	Width	Length	Rolls/Case	
3M™ Scotchcast™ Plus Casting Tape					
82001	White	1" (2,5 cm)	2 yds. (1,8 m)	10	
8200*	White	2", 3", 4", 5" (5 cm, 7,5 cm, 10 cm, 12,5 cm)	4 yds. (3,6 m)	10	
8200*A	Black	2", 3", 4" (5 cm, 7,5 cm, 10 cm)	4 yds. (3,6 m)	10	
8200*B	Blue	2", 3", 4", 5" (5 cm, 7,5 cm, 10 cm, 12,5 cm)	4 yds. (3,6 m)	10	
8200*G	Green	2", 3", 4", 5" (5 cm, 7,5 cm, 10 cm, 12,5 cm)	4 yds. (3,6 m)	10	
8200*J	Glow-in-the-dark	2", 3", 4" (5 cm, 7,5 cm, 10 cm)	4 yds. (3,6 m)	10	
8200*L	Light Blue	2", 3", 4", 5" (5 cm, 7,5 cm, 10 cm, 12,5 cm)	4 yds. (3,6 m)	10	
8200*R	Red	2", 3", 4", 5" (5 cm, 7,5 cm, 10 cm)	4 yds. (3,6 m)	10	
8200*U	Purple	2", 3", 4", 5" (5 cm, 7,5 cm, 10 cm)	4 yds. (3,6 m)	10	
8200*V	Bright Green	2", 3", 4" (5 cm, 7,5 cm, 10 cm)	4 yds. (3,6 m)	10	
8200*W	Bright Orange	2", 3", 4" (5 cm, 7,5 cm, 10 cm)	4 yds. (3,6 m)	10	
8200*X	Bright Pink	2", 3", 4" (5 cm, 7,5 cm, 10 cm)	4 yds. (3,6 m)	10	
8200*Y	Yellow	2", 3" (5 cm, 7,5 cm)	4 yds. (3,6 m)	10	
8200*M	Maroon	2", 3" (5 cm, 7,5 cm)	4 yds. (3,6 m)	10	
3M™ Scotchcast™ Soft Cast Casting Tape					
82101	White	1" (2,5 cm)	2 yds. (1,8 m)	10	
8210*	White	2", 3", 4", 5" (5 cm, 7,5 cm, 10 cm, 12,5 cm)	4 yds. (3,6 m)	10	
8210*B	Blue	2", 3", 4" (5 cm, 7,5 cm, 10 cm)	4 yds. (3,6 m)	10	
8210*U	Purple	2", 3", 4" (5 cm, 7,5 cm, 10 cm)	4 yds. (3,6 m)	10	
8210*R	Red	2", 3", 4" (5 cm, 7,5 cm, 10 cm)	4 yds. (3,6 m)	10	
3M™ Scotchcast™ Wet or Dry Cast Padding				Each/Bag	Bags/Case
WDP2		2" x 4 yds. (5 cm x 3,6 m)	20	4	
WDP3		3" x 4 yds. (7,5 cm x 3,6 m)	20	4	
WDP4		4" x 4 yds. (10 cm x 3,6 m)	20	4	
WDP6		6" x 4 yds. (15 cm x 3,6 m)	20	4	

Continued on back page

3M™ Scotchcast™ Casting Tapes

Celebrating over 30 years of innovation



Innovation in Immobilization

3M

3M and Scotchcast are trademarks of 3M.
Please recycle. Printed in U.S.A.
© 3M 2014. All rights reserved.
70-2009-4217-8

3M

3M™ Scotchcast™ Family of Casting Tapes

3M™ Scotchcast™ Casting Tapes allow you to give your patients the mobility they need to return to most normal daily activities, as their injury allows. Our family of tapes provides the choice of fiberglass or polyester backing. Our casting products are supported by a full line of accessories.

3M™ Scotchcast™ Plus Casting Tape

A lightweight, strong and durable casting tape that combines the benefits of a fiberglass casting tape with the handling ease of plaster. Available in 13 colors to meet patient preferences.



Product Benefits:

- Water-activated resin lets you control working time according to application needs.
- High and early strength allows for weight bearing in 20 minutes and provides a durable cast with less chance of breakdown and the need for costly repairs or replacements.
- Lightweight for easy mobility.
- Stretches to conform well to contours for a more comfortable, better fitting cast.
- Porous weave allows the skin to breathe.
- Available in 13 colors to maximize patient satisfaction.
- Slippery resin formula eliminates the need for cast creams or special gloves.

Suggested Applications:

- For use in the treatment of musculoskeletal injuries.

Available Colors:

- | | | |
|--------------------|------------------------|-------------------|
| ■ Black (A) | □ White | ■ Blue (B) |
| ■ Green (G) | ■ Light Blue (L) | ■ Red (R) |
| ■ Purple (U) | ■ Bright Orange (W) | ■ Bright Pink (X) |
| ■ Bright Green (V) | ■ Glow in the Dark (J) | ■ Yellow (Y) |
| ■ Maroon (M) | | |

3M™ Scotchcast™ Wet or Dry Cast Padding

For use with 3M™ Synthetic Cast Stockinet and 3M™ Scotchcast™ Plus Casting Tape. A hypoallergenic, synthetic cast padding that allows patients to get their cast wet, when following the instructions for use.



Product Benefits:

- Lets you allow patients to shower, bathe, swim or use a whirlpool.
- Sheds water and feels comfortably dry to patient in approximately 1 to 3 hours.*
- Allows you to stock only one cast padding for all your patient needs.
- Suitable for use with all synthetic casting materials such as 3M™ Scotchcast™ Plus casting tape.

Suggested Applications:

- For use in the treatment of musculoskeletal injuries.



3M™ Scotchcast™ Soft Cast Casting Tape

A fiberglass casting tape that does not set to a completely rigid state, but remains flexible and soft even when it is completely cured.



Product Benefits:

- Flexible, soft cast allows for mobility.
- The soft, pliable cast edges improve patient comfort.
- Can be used in combination with rigid fiberglass splints, as soft edges provide additional comfort.

Suggested Applications:

- For use in the construction of casts where support is needed but rigid immobilization is not required, or a limited range in motion is desired.

Applications include:

- Sports Medicine: thumb, wrist and ankle sprains
- Pediatric Orthopedics: serial casting for club foot treatment
- General Orthopedics: secondary casting, functional stabilization casting
- Occupational Therapy: removable splints
- Diabetic Foot: total contact casting

Available Colors:

- | | |
|--------------|------------|
| □ White | ■ Blue (B) |
| ■ Purple (u) | ■ Red (R) |

3M™ Casting Accessories



3M™ Scotchcast™ Wet or Dry Cast Padding

A hypoallergenic, synthetic cast padding that allows patients to get their cast wet when following the instructions for use. Suitable for use with all synthetic casting materials, such as 3M™ Scotchcast™ Plus casting tape.



3M™ Synthetic Cast Padding

A non-woven polyester padding that is porous, highly conformable, strong, yet easy to tear. Suitable for use with all casting materials, plaster and synthetic.



3M™ Synthetic Cast Stockinet

Essential component of a high-quality cast, this synthetic stockinet is intended for use as an underlayment for all standard casting applications. Suitable for use with all casting materials, plaster and synthetic.


Application Instructions for Scotchcast™ Plus and Poly Tapes



Apply stockinet and cast padding (gloves must be worn during application of casting tapes).



Remove roll from foil pouch and dip in cool water, squeezing once. With exposure to more water, the setting time will be faster than with exposure to less water.

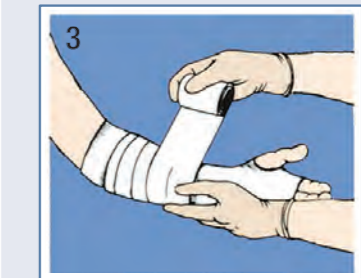

Apply in spiral motion, over-lapping by half on each turn.



Smooth and mold cast to shape of extremity. Weight bearing in 20 minutes.

Call your 3M Health Care Representative for specific uses and application techniques for 3M™ Scotchcast™ Soft Cast Casting Tapes.