


 **smith&nephew**
MELOLITE[◇]
Low Adherent Wound Dressing

Reduces trauma to healing tissue

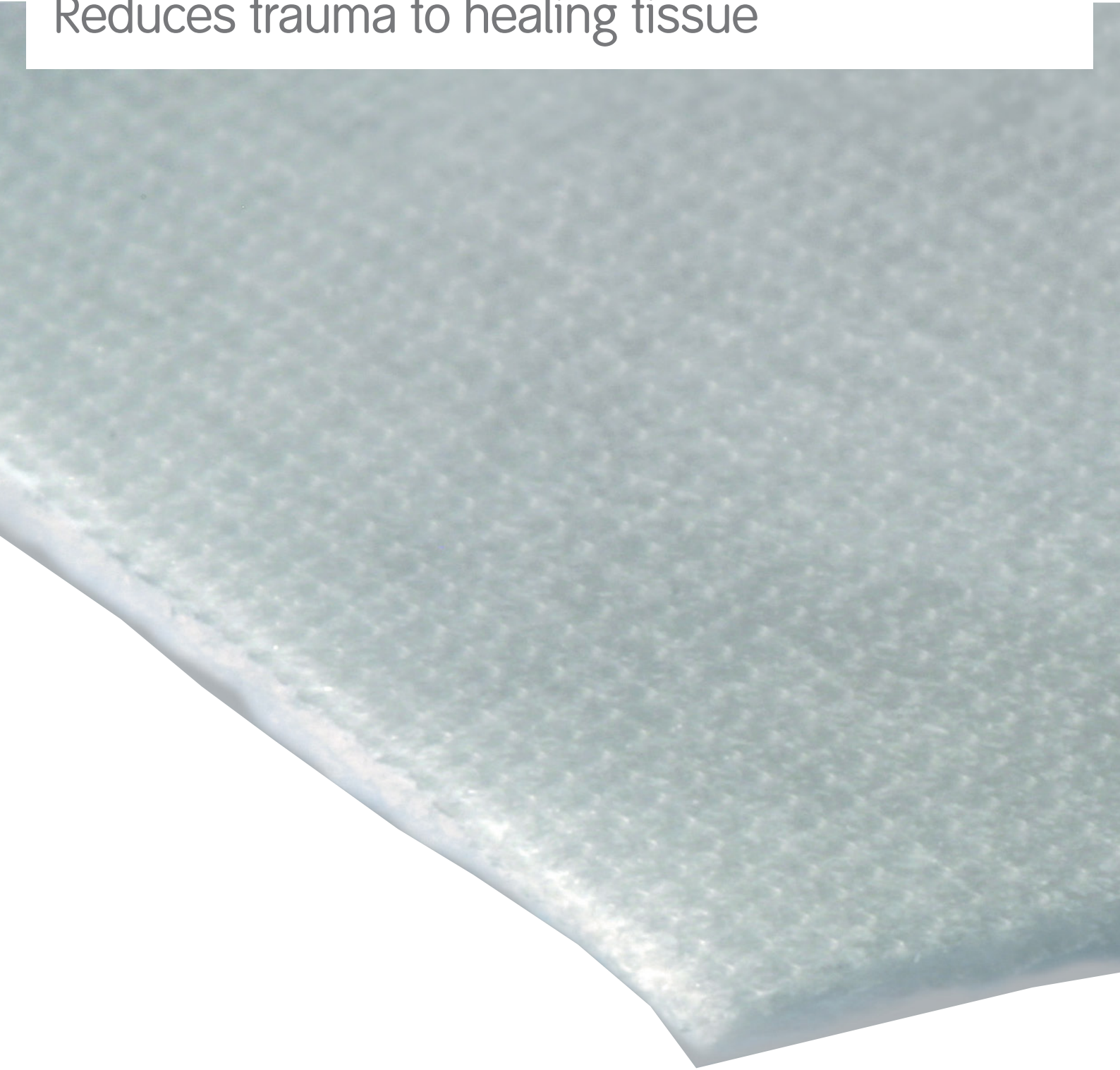

MELOLITE[◊]: Low adherent absorbent dressing

MELOLITE[◊] dressings are constructed from an absorbent pad which helps distribute the fluid, thus avoiding the problems of immediate exudate strike-through.

Advantages

Minimises risk of wound trauma on dressing removal

MELOLITE is covered on both sides with a heat-laminated polyethylene film which acts as a low-adherent wound contact layer by allowing the passage of exudate through its elliptical holed construction into the absorbent pad. This minimises the risk of adherence to the wound and the related trauma on dressing removal.

Easy to use

Either side of the dressing can be placed directly onto the wound site reducing time spent in dressing changes.

Patient comfort

The softness of the dressing provides patient comfort¹ and high conformability to wound sites.

Indications


MELOLITE is indicated for casualty and postoperative use over clean, sutured wounds and for lacerations and abrasions

Presentation

Each MELOLITE dressing is individually wrapped and sterile.

Ordering information

Code	Product description (dressing; pad)	Items per unit
4811	7.5cm x 5cm	Box/100
4812	7.5cm x 10cm	Box/100
4814	7.5cm x 20cm	Box/100

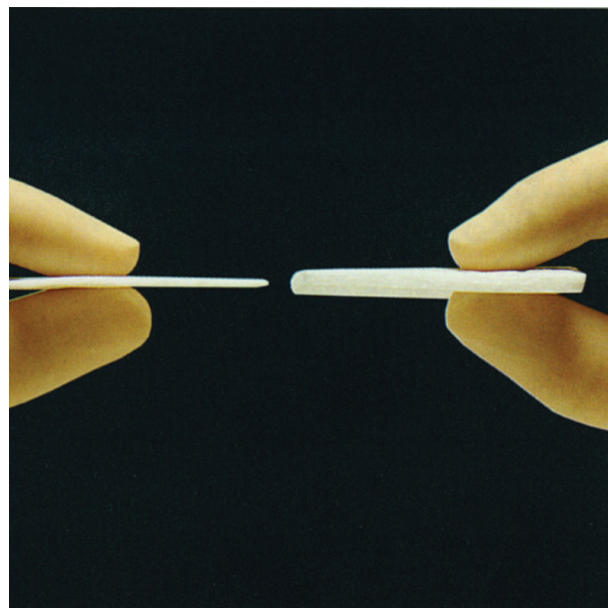

References

1. Parish HA, Roberts CD. An evaluation of a lightweight, non-adherent dry dressing in the clinical situation. Data on file; 1983: protocol UK 34.

Smith & Nephew Pty Ltd Australia

T +61 2 9857 3999 F +61 2 9857 3900
www.smith-nephew.com/australia

Customer Service

T 13 13 60 F 1800 671 000

Smith & Nephew Ltd New Zealand

T +64 9 820 2840 F +64 9 820 2841
www.smith-nephew.com/new-zealand

Customer Service

T 0800 657 799 (Surgical) T 0800 807 663 (Wound)

Please consult product labels and inserts for any indications, contraindications, hazards, warnings, cautions and instructions for use.

[◊]Trademark of Smith & Nephew
SNI2763 (11/16)